

The Promise

White Laser Healing System

with Tachyonized Pyramidal Photoporation Chamber



White Laser

Instruction Manual

Setting up instruction

For amplified Crystal therapy:

1. Put the fiber optic cable to the laser aperture.



2. Remove copper head by twisting and pulling out from the wand.



3. Place crystal bead of choice into the pyramid head and place back the copper head.



4. Seven beads resonating with the 7 chakras and their associated organs.

See color / region chart in treatment manual

5. Clear crystal use dropper or syringe to place drops of remedy into clear crystal and seal with silicon cap tightly.



6. Place remedy filled clear crystal into the freezer.

7. After freezing, place remedy filled clear crystal into head as above and start the projection pattern and expose to required point according to the part being healed. See point's diagram in separate manual.

For keeping geometry shapes energy:

1. Remove the fiber optic cable.



2. Remove copper head by twisting and pulling out from the wand.



3. Place crystal bead of choice into the pyramid head and place back the copper head.



4. Put the copper crystal head into the laser aperture.



5. Clear crystal use dropper or syringe to place drop of remedy into clear crystal and seal with silicon cap tightly.

6. Place remedy filled clear crystal into head as above.

7. Remove copper head again and take the crystal out was and store in box keep it with you.



Attention: 1. Please, place the laser projector in the bare space.

2. Please, unplug the laser projector when it is not used for a long time.

Project the flower remedies mood through the crystal.

To create personal programmed crystals place the crystal to be programmed in the holder and start the pattern required to begin the programming. Remove the crystal and put it in silver holder around your neck.

To erase crystals you require a Tachyon Photon Wand to clear the crystal. If you record with light it's stored there and you can carry with you. If you take 2 crystals cut from the same rock one on the projector the other on your neck, you receive full power of projection while away from your place. Can also use in remote healing by sending the programmed crystal to the remote person to establish the connection.

Operation Instruction

1. Plug in the electric wire in the back of the laser projector.



2. Put the SD card which comes with the White Laser in the back of the laser projector.



3. Turn the switch “ON or OFF” to turn on or turn off.



4. Using the key to switch on or shut down the laser irradiation.



Attention: Please “shut down” the laser irradiation before turn off the laser projector.


5. Choose the folder of the geometry shape of choice by twisting the shuttle button and press.

Attention: Choose the icon  going back to the folder list.

6. Choose the file of your choice.

Attention: The shape can be expand and drop by twisting the small button.

Main Interface Function Description

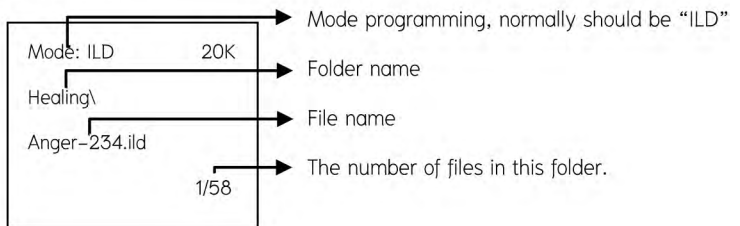

Contraindications

1. Pregnant women
2. Cancer
3. Hemorrhagic diseases

Troble–shooting instructions

| Faults | Causes | Solutions |
|--|------------------------------------|--|
| Cannot turn on | No electric power | Check the electric wire. |
| To select English Language | | <ol style="list-style-type: none">1. Turn off the projector2. Turn on the projector again while the loading screen appear press the shuttle button and hold for 2 seconds.3. Select “English” and press the shuttle button again |
| Projector on, but no laser output | | Insert and turn key to the on position. |
| Program menu not displayed | SD card not fully inserted. | Press the SD card into the host tightly. |
| | Problem with SD card loading. | Power the projector off then on again. Attention: If still unsolved, please contact our tech support. |
| Message on screen asks for list file and doesn't load the SD card menu | The projector play mode was reset. | <ol style="list-style-type: none">1. Push the shuttle button for 3 seconds.2. Choose the “Mode select” and Choose ILD file. |

Maintenance & cleaning

1. Use cotton soaked in neutral solution to clean the surface. Do not use organic solution or gasoline.
2. Please using cotton balls soaked in alcohol to clean the optic part and probe. Do not clean it with something abrasive.
3. Please keep away from humidity, heat source and direct exposure to sunlight, in order to avoid the degradation of components which may shorten the service life.

Product warranty

Dear Customers:

Thank you for your purchase of The Promise Laser Product, your satisfaction is our priority. Please register here once you receive the product to activate the warranty. The Promise, will provide free service warranty for one year against manufacture defects.

guarantees;

1. Damages caused by improper operation or mistakes.
2. Breakdowns caused by self-disassembling.
3. Material costs should be charged for products over the warranty period.
4. Breakdowns or damages caused by accidents or human activities (dropping, scratches,

soaking or self-replacement of parts).

Product List

| | | |
|---------------------------|---|---------------|
| 1. Laser projector | 1 | piece |
| 2. Laser optical wire | 1 | piece |
| 5. Crystal wand | 1 | piece |
| 6. Crystal set | 1 | box (8 beads) |
| 7. Operation instructions | 1 | book |
| 8. Function instructions | 1 | book |

Laser projector specification

| | |
|------------------------------|----------------------|
| Model | 850mW |
| Dimensions and Weight | |
| Length x Width x Height | 316 x 254 x 150 mm |
| Weight | 4.26 Kg |
| Power | |
| Power Supply Type | Switch in (Internal) |
| Range | 100–240V 50/60Hz |
| Consumption | 26W |
| Fuse | T1.6A 250V |
| Light source | |
| Red | 638nm 200mw |
| Green | 532nm 150mw |
| Blue | 450nm 500mw |
| Type | DPSS Class 4 Laser |
| Scanner speed | 20K |
| Working condition | |
| Maximum external Temp | 40 °C |
| Cooling System | Green TEC cooling |